



easyCS™ Circular Stapler




- **Smaller head diameter**
Smaller head diameter of the anvil provides easier control for insertion and removal during placement.
- **Shorter anvil height**
Shorter anvil height design benefits the kinking anatomy of esophagus for better solution.

- **Shorter staple housing height**
Shorter staple housing height is designed for increased ability of access to the narrow space and decreased incidence of injury to intestinal wall.



Product Information

Code	Open Staple Height	Color	Anvil Diameter
CS2135	3.5mm	Sky blue	21mm
CS2148	4.8mm		
CS2535	3.5mm	White	25mm
CS2548	4.8mm		
CS2935	3.5mm	Dark blue	29mm
CS2948	4.8mm		
CS3148	4.8mm	Green	31mm
CS3348	4.8mm	Yellow	33mm

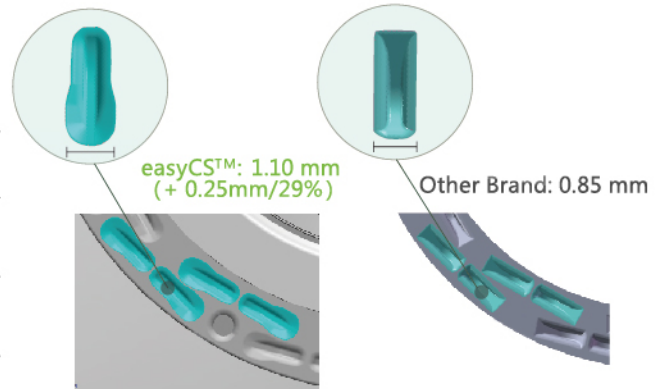

Tissue contacting surface will greatly guide staples into staple pocket which enables leak-resistant and hemostatic anastomosis.

EZISURG MEDICAL CO., LTD.

Add: Bldg.10, 588 Tianxiong Rd., SIMZ, Pudong District, Shanghai, China
 Tel : +86-21-5045 6176 Fax : +86-21-5067 6156
 Email : info@ezisurg.com
 www.ezisurg.com



CE 0197 EN ISO 13485:2012