
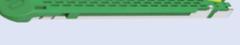


## Product Information

	Code	Shaft Length	Shaft OD	Package(BX)
	U12S45	260mm	12.3mm	1
	U12M45	350mm	12.3mm	1
	U12L45	440mm	12.3mm	1
	U12S60	260mm	12.3mm	1
	U12M60	350mm	12.3mm	1
	U12L60	440mm	12.3mm	1

Product Series	Product Code	Cartridge Color	Number of Staple Rows	Open Staple Height (mm)	Closed Staple Height (mm)	EA/BX
	N45M	Grey	6	2.0	0.75	6
	N45W	White	6	2.5	1.0	6
	N45B	Blue	6	3.5	1.5	6
	N45C	Gold	6	3.8	1.7	6
	N45G	Green	6	4.1	2.0	6
	N60M	Grey	6	2.0	0.75	6
	N60W	White	6	2.5	1.0	6
	N60B	Blue	6	3.5	1.5	6
	N60C	Gold	6	3.8	1.7	6
	N60G	Green	6	4.1	2.0	6
	N45BC	Blue	6	3.5	1.5	6
	N60GC	Green	6	4.1	2.0	6
	N60CC	Gold	6	3.8	1.7	6
	N60BC	Blue	6	3.5	1.5	6
	N45GC	Green	6	4.1	2.0	6
	N45CC	Gold	6	3.8	1.7	6
	N45VM	Brown	6	3.5-2.5-2	1.5-0.88	6
	N45MT	Purple	6	4.1-3.8-3.5	2.25-1.5	6
	N45XT	Black	6	4.4-4.1-3.8	3.0-2.25	6
	N60VM	Brown	6	3.5-2.5-2	1.5-0.88	6
	N60MT	Purple	6	4.1-3.8-3.5	2.25-1.5	6
	N60XT	Black	6	4.4-4.1-3.8	3.0-2.25	6

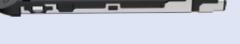

### EZISURG MEDICAL CO., LTD.

Add: No.399, Miaoqiao Rd., Pudong New District, 201315 Shanghai, P.R. China  
 Tel: +86-21-5045 6176 Fax: +86-21-5067 6156  
 Email: [contact@ezisurg.com](mailto:contact@ezisurg.com)  
 Website: [en.ezisurg.com](http://en.ezisurg.com)




# easyEndo™ Lite

## Linear Cutting Staplers and Loading Units for Single Use

## 120 degrees of articulation


A smoother and more precise positioning in various clinical scenarios, improving efficiency during operation.

**easyEndo™ Lite allows the surgeon to rotate and articulate in 9 different angles.**

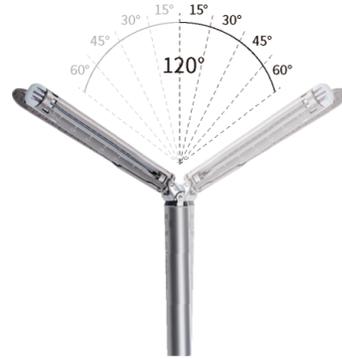

### Greater Mobility

## easyEndo™ Lite 120° articulation in bariatric surgeries

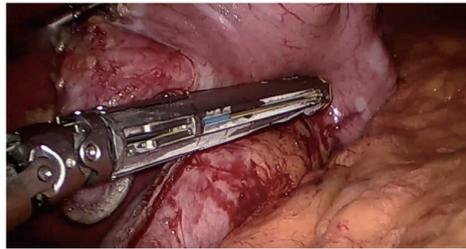

Easier to perform vertical stomach cuts in Sleeve surgery, ROUX-EN-YGASTRIC BYPASS, and in crowded spaces such as OAGB.


## easyEndo™ Lite Sharpest knife



easyEndo™ Lite has improved the hardness of the knife blade (**500Hv**) to reduce hemorrhage risk.

Industry standard (377)

## Smart Cartridge

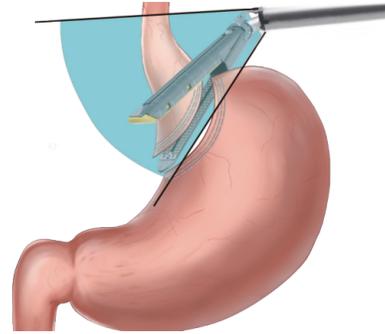

### Thinner fixed anvil

Thinner Fixed Anvil provides an increased level of agility and allows the user to operate in small and narrow spaces inside the body.

## easyEndo™ Lite has reduced the size of the Cartridge slot




## Ergonomic design

easyEndo™ Lite double release mechanism reduces the strain on the wrist and hands of surgeons.




## Better clamping force than previous easyEndo™ Linear Cutters

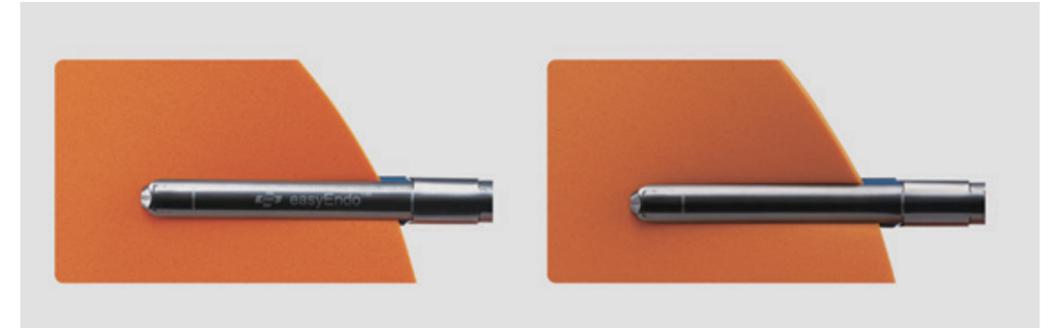

easyEndo™ Lite has reduced tissue spillage during the firing process, preventing tissue injury caused by over pressing.

## Intelligent Triangular-Anvil Shaped Bump

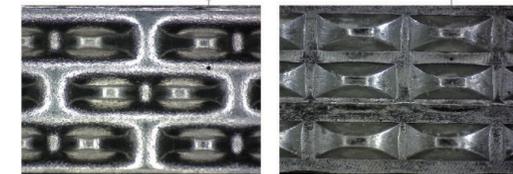

ITASB - Technology


No ITASB - Technology

Fully effective staple guiding and compound anvil technology to ensure that staples have a better B-shaped formation.

100% of the receiving area for staple leg guiding provides a perfect B-shaped formation.

Perfect B-shaped formation to achieve the best suture and hemostasis effect.

This speeds up the healing process by providing better blood perfusion.

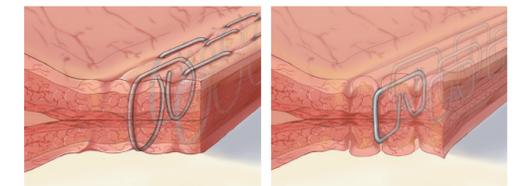

B - Shaped Formation