
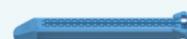
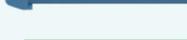
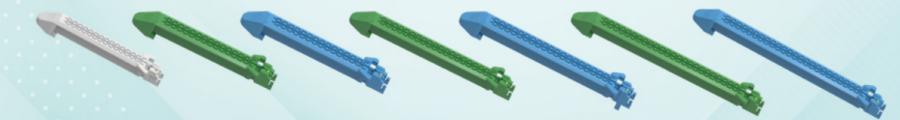


3+3 Model Options

Stapler	Model	Pre-loading Unit	Color of Cartridge	Staple Line Length(mm)
	LCE60X	/	/	60
	LCE60	LCC60W	● White	60
	LCE60B	LCC60B	● Blue	60
	LCE80X	/	/	80
	LCE80	LCC80B	● Blue	80
	LCE80G	LCC80G	● Green	80
	LCE100B	LCC100B	● Blue	100
	LCE100	LCC100G	● Green	100

Cartridge	Model	Suture Length (mm)	Stapler Height (mm)	Closed Height (mm)	Color
	LCC60W	60	2.5	1.0	● White
	LCC60B	60	3.8	1.5	● Blue
	LCC60G	60	4.8	2.0	● Green
	LCC80B	80	3.8	1.5	● Blue
	LCC80G	80	4.8	2.0	● Green
	LCC100B	100	3.8	1.5	● Blue
	LCC100G	100	4.8	2.0	● Green

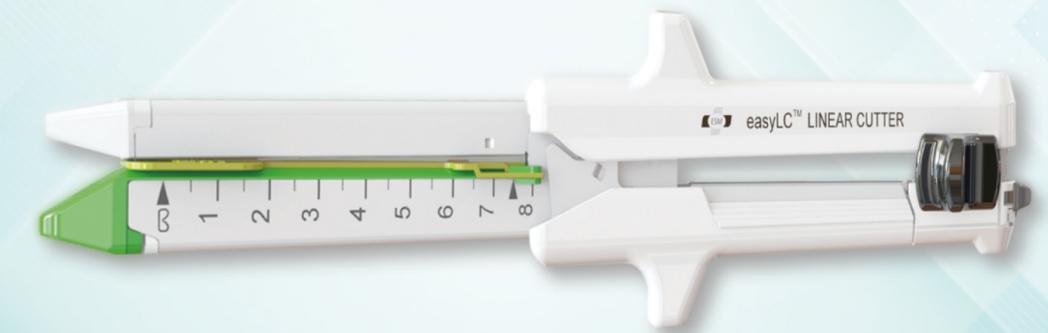

easyLC™

Linear Cutting Staplers
and Loading Units for Single Use

Procedures:



- Colectomy
- Small Bowel Resection
- Gastrectomy
- Lobectomy Wedge Resection

Specialties:



- Gastrointestinal
- Gynecologic
- Thoracic

Optimized B-shape formation


With the **optimized B-shape formation technology**, this arc-shaped staple ensures more secure compression and greater blood flow, promoting tissue healing.

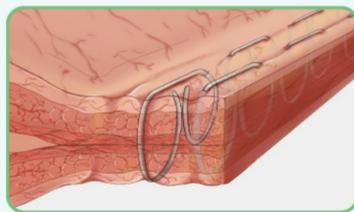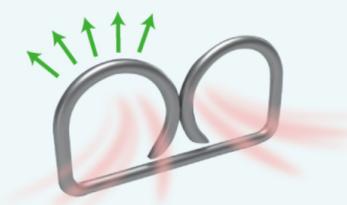

easyLC™

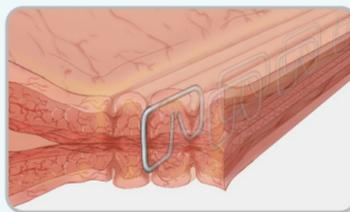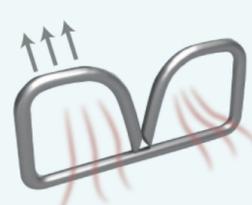


Others

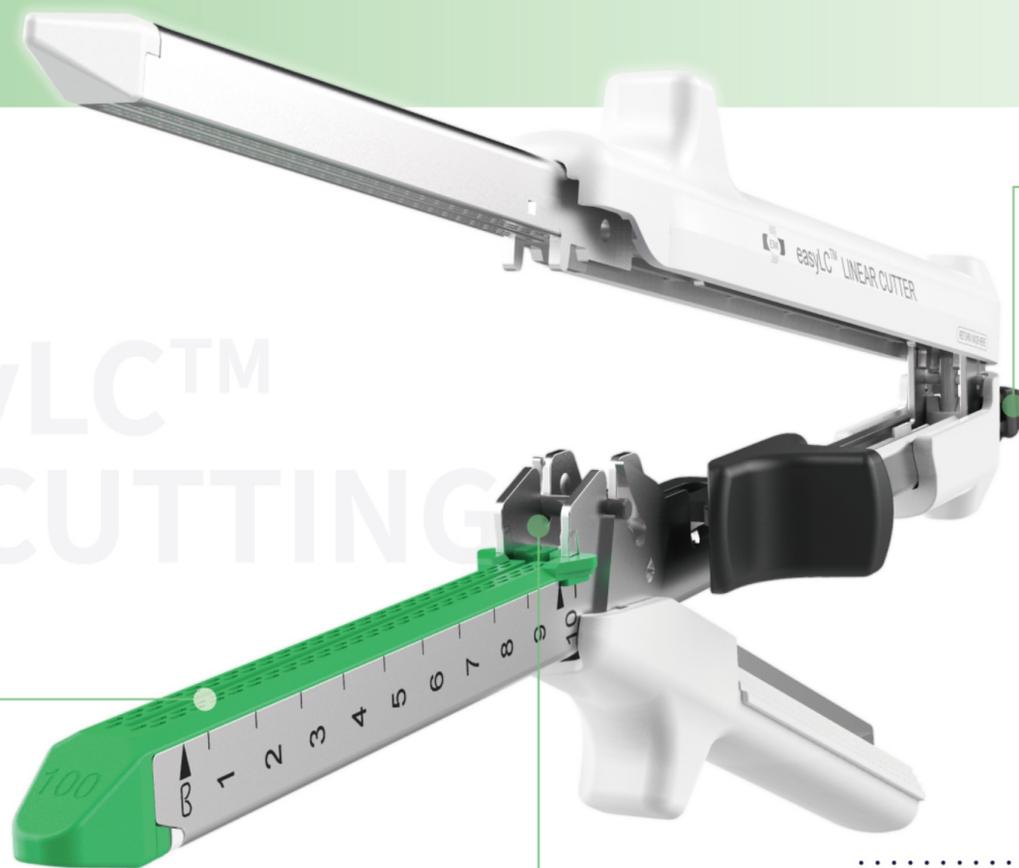


easyLC™
LINEAR CUTTING

Tempered stainless steel blade

The easyLC™'s stainless steel blade is at least **20% harder** than comparable products, which helps achieve smoother and more precise cutting, reducing bleeding caused by tearing.

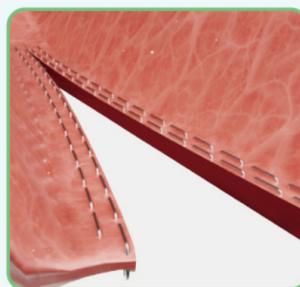

easyLC™ (blade rigidity: 500Hv)

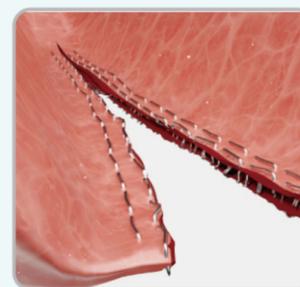


Others (blade rigidity: ≤400Hv)


Adjustable clamping force

The easyLC™ offers **3 levels of clamping force options**, making it suitable for different tissue types. This feature aids in hemostasis and helps prevent tissue necrosis.

1 Level 1 for thinner tissue

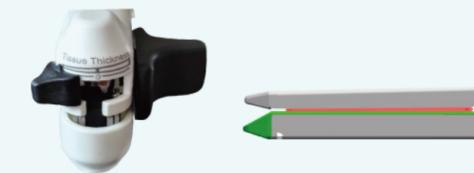

2 Level 2 for common tissue

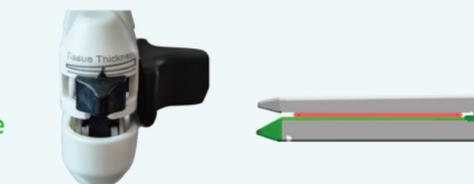

3 Level 3 for thicker tissue

